

Anabolic Lab™

Analytical Lab Testing for
Anabolic Steroid Harm Reduction

<https://anaboliclab.com>

Lab Report – May 4, 2016

Med-Tech Solutions
Stanabol (10 mg/tablet)
Stanozolol



If you like the attached report, please visit [our website](https://anaboliclab.com) to find out how you can support independent and objective lab testing.



Simec AG

Areal Bleiche West Postfach 413 CH - 4800 Zofingen
Telefon +41 62 752 83 08 Fax +41 62 752 83 09
E-Mail: info@simec.ch www.simec.ch

Test report

Number **C 42373**

Object of analysis: Med-Tech Stanabol

Customer: AnabolicLab.com

Samples: 1

Batch: 0001

Receipt: 24.03.2016


Subject: Assay, Identification Stanozolol

Method(s): HPLC

Sample / Analysis:	Method	Status	Result	Specification	complies
<u>Stanozolol:</u>					
Identification	HPLC	N	complies	complies	yes
Content	HPLC	N	9.53 mg/tablet	10 mg/tablet	-

Method-status: G = GMP A = accredited V = generally validated P = validated on product N = not validated
E = external lab

Remarks: The sample(s) mentioned above have been analysed as they have been sent to us by the client. We have not controlled external processes like manufacturing, labelling, sampling, shipping, storage. These results are for information only and do not compensate for correct quality control by the manufacturer/distributor. Generally, pharmaceutical products have to be produced and distributed under full GMP/GDP regime, including GMP-compliant analyses with validated methods.

Signatures: controlled & released: 

Staff: 

completed:
04.05.2016


The enclosed results refer exclusively to the object of examination described above. The measuring accuracy is available on request. Only handwritten signatures are valid. Without written permission, it is not allowed to publish single parts of this report.

page 1 of 1

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 04.05.2016 06:15:40
 Vial: 90
 Probenidentifikation: C42373 A
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 121, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 04.05.2016 11:38:30
 Injektvolumen: 20

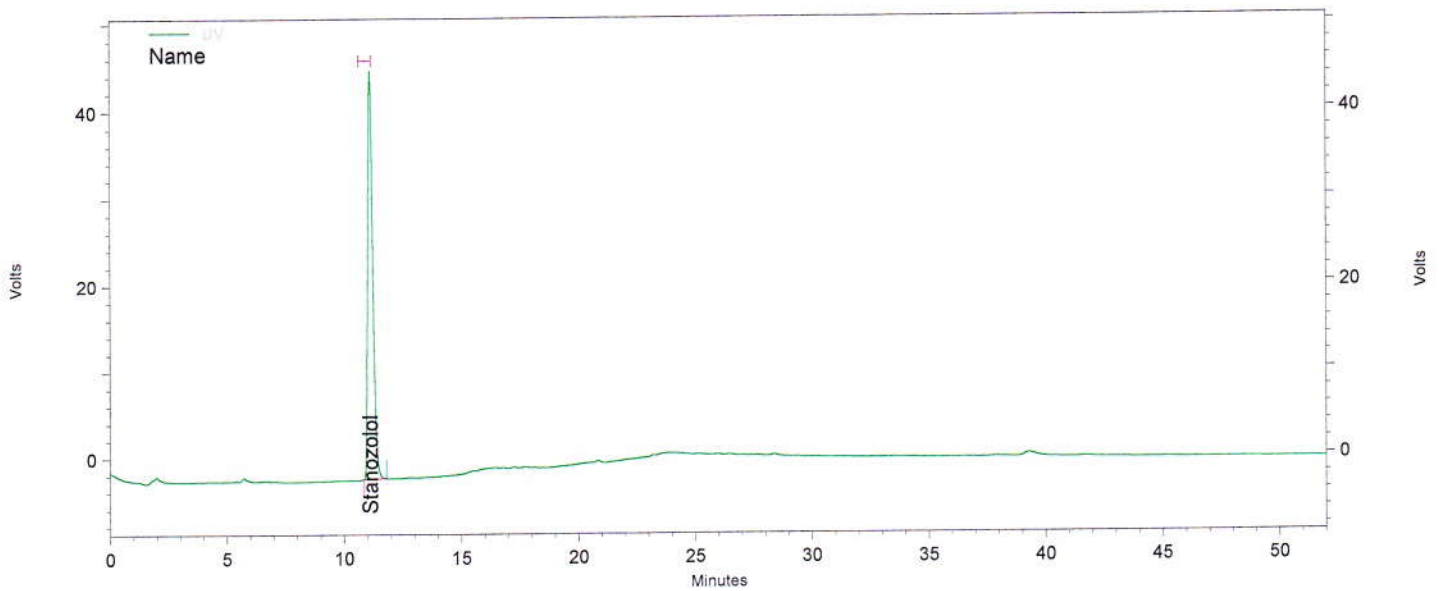

Channel A


UV Results				
Name	Time	Area	Amount	Area Percent
Stanozolol	11,120	2790351	0,000	100,000
Methenolone enanthate			0,000 BDL	
Totals		2790351	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse
 Steroide\019_C42373 A-Rep1.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 04.05.2016 07:09:06
 Vial: 90
 Probenidentifikation: C42373 A
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 121, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 04.05.2016 11:38:34
 Injektvolumen: 20

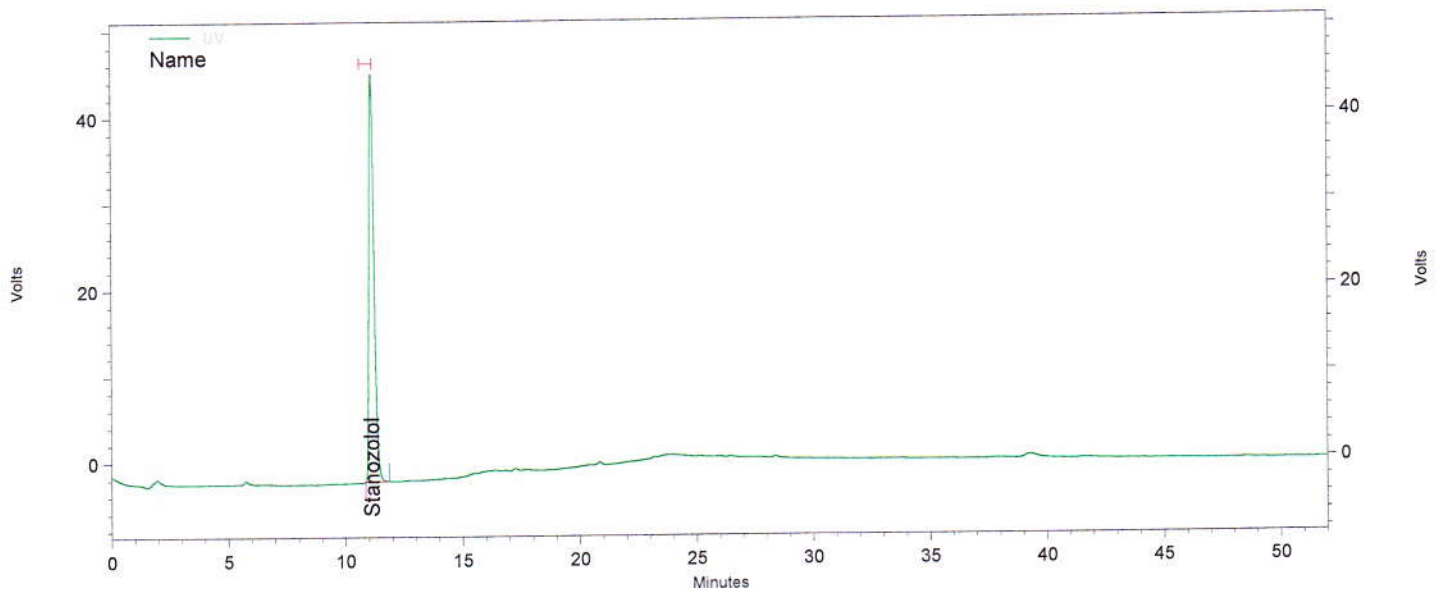

Channel A


UV Results				
Name	Time	Area	Amount	Area Percent
Stanozolol	11,117	2800352	0,000	100,000
Methenolone enanthate			0,000 BDL	
Totals		2800352	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse Steroide\019_C42373 A-Rep2.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 04.05.2016 08:02:32
 Vial: 90
 Probenidentifikation: C42373 A
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 121, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 04.05.2016 11:38:38
 Injektvolumen: 20

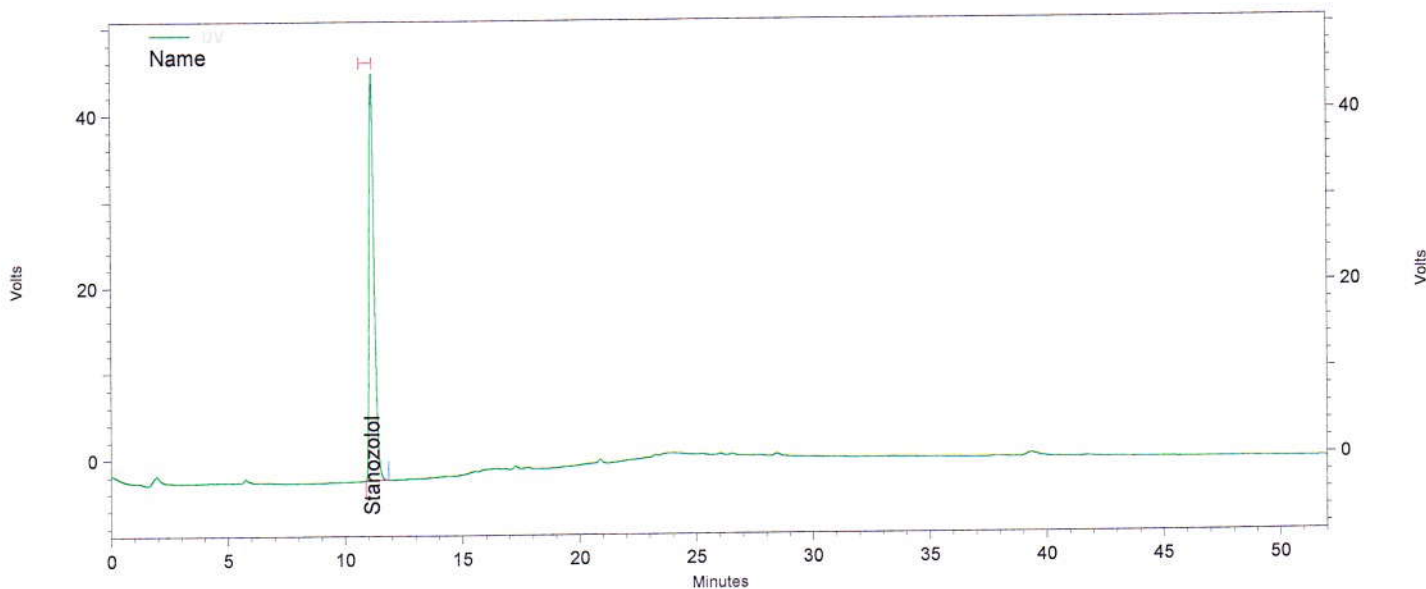

Channel A


UV Results				
Name	Time	Area	Amount	Area Percent
Stanozolol	11,143	2781178	0,000	100,000
Methenolone enanthate			0,000 BDL	
Totals		2781178	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse
 Steroide\019_C42373 A-Rep3.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 04.05.2016 08:55:56
 Vial: 91
 Probenidentifikation: C42373 B
 Operator: admin
 Chrom.-Beding.: ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule: Nr. 121, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 04.05.2016 11:38:42
 Injektvolumen: 20

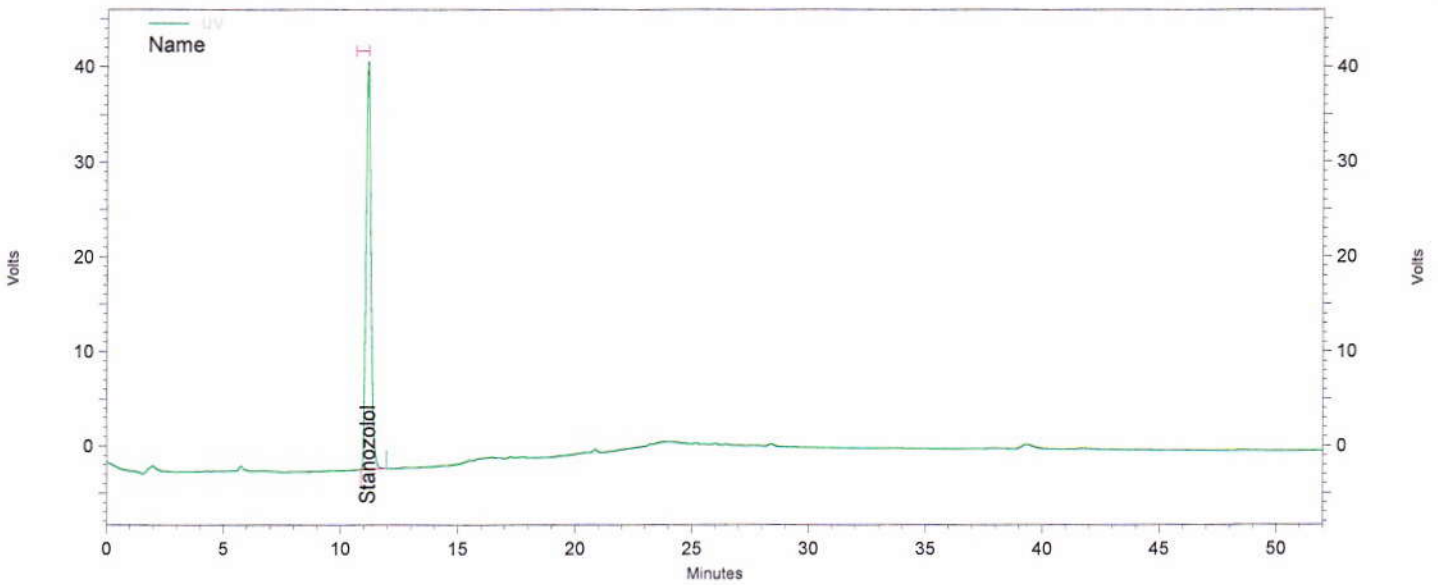

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Stanozolol	11,140	2513586	0,000	100,000
Methenolone enanthate			0,000 BDL	


Totals		2513586	0,000	100,000
--------	--	---------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse Steroide\020_C42373 B-Rep1.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 04.05.2016 09:49:20
 Vial: 91
 Probenidentifikation: C42373 B
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 121, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 04.05.2016 11:38:46
 Injektvolumen: 20

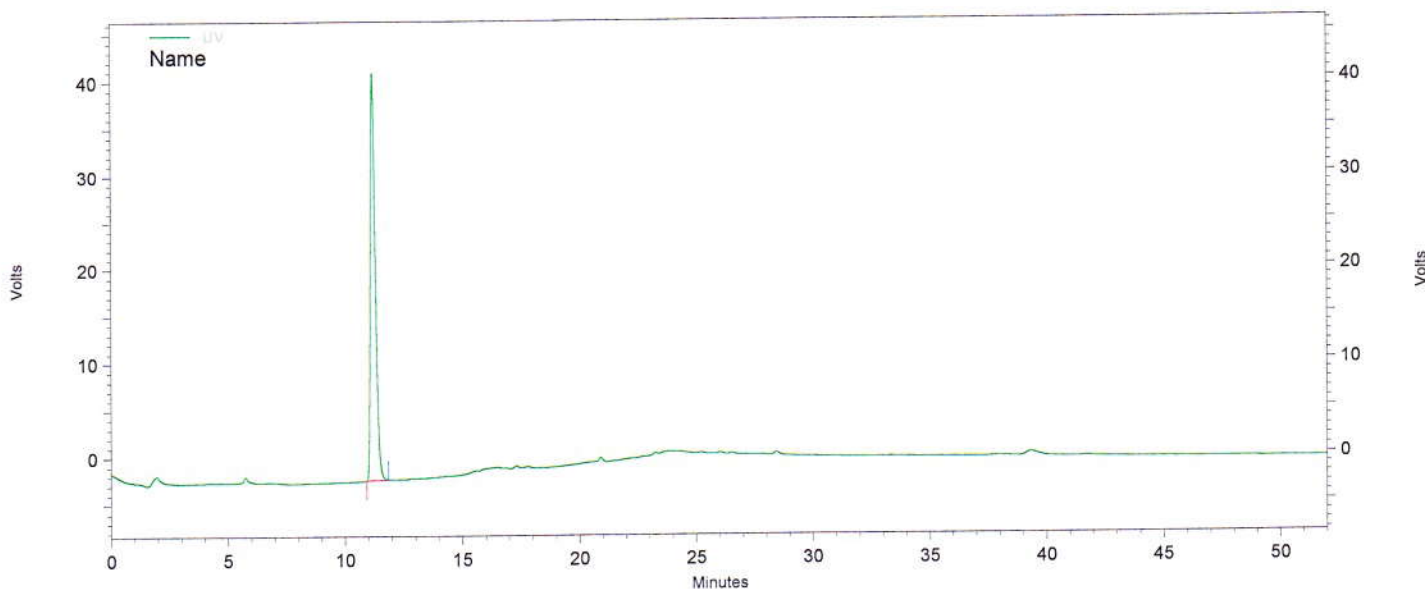


Channel A

UV Results				
Name	Time	Area	Amount	Area Percent
Stanozolol	11,187	2493714	0,000 BDL	100,000
Methenolone enanthate			0,000 BDL	
Totals		2493714	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse Steroide\020_C42373 B-Rep2.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 04.05.2016 10:42:42 Druckdatum: 04.05.2016 11:38:51
 Vial: 91 Injektvolumen: 20
 Probenidentifikation: C42373 B
 Operator: admin
 Chrom.-Beding.: ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule: Nr. 121, 250 x 4,6, C18
 Detektion: UV 240nm

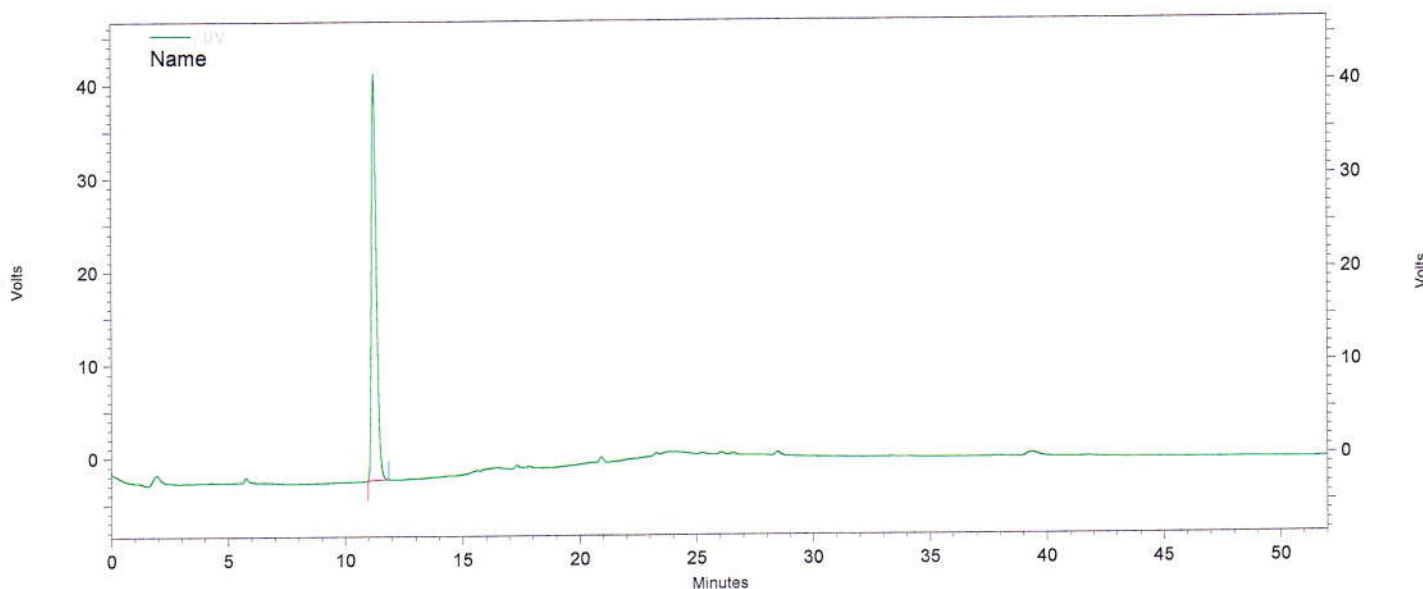

Channel A

UV Results				
Name	Time	Area	Amount	Area Percent
Stanozolol	11,217	2486225	0,000 BDL	100,000
Methenolone enanthate			0,000 BDL	
Totals		2486225	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\Mai 2016\02,05,16\Diverse Steroide\020_C42373 B-Rep3.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met